

## BURNERS, MEKER

Nickel plated brass tube fitted with deep ceramic grid 24mm diameter to give a high and evenly distributed flame temperature. Stove enamel steel base.

	Gas type
<b>C21-2700</b>	Natural
<b>C21-2750</b>	LPG

## BURNERS, MEKER

As **C21-2700** series but with 30mm diameter ceramic grid.

	Gas type
<b>C21-3000</b>	Natural
<b>C21-3050</b>	LPG

## BURNERS, BUNSEN, UNIVERSAL

Nickel plated brass tube, 13mm diameter, air regulator, needle valve control and adjustable baffle stabiliser. For natural gas, LPG, coal gas or mixtures of these gases.

**C21-3200** Universal Bunsen Burner

## FLAME SPREADER

Nickel plated brass, 50mm wide at top.

	To fit burners
<b>C21-3300</b>	Tube dia. 11mm, natural gas and LPG
<b>C21-3310</b>	Tube dia. 13mm, natural gas and LPG

## BURNER, RING

Cast iron burner ring mounted on a heat deflecting vitreous enamelled base, with air and gas regulators, suitable for natural gas only. Dimensions 210 x 300 x 70mm (dia. x l x h).

**C21-3400** Gas Burner

## BURNER, CARTRIDGE

Designed for use with Butane/Propane mix gas cartridges (NOT supplied). Brass construction with replaceable plated steel perforated head.

**C21-3500** Cartridge Burner

## BURNER, SPIRIT

Glass spirit lamp with fitted wick in ceramic burner, nominal capacity 100ml. Pack of 10

<b>C21-3610</b>	Spirit Burner
<b>C21-3650</b>	Wick, spare, 10m

## BURNER, SAFETY

Electric ignition, activated by hand sensor or foot switch accessory. For use with natural gas, LPG or coal gas; gas/air mixture adjustable. Automatic gas cut-off with over-heating, non-ignition or failure. Complete with built-in timer and 3 interchangeable gas nozzles. For 110/220V a.c. single phase mains supplies and connection to a permanent gas supply.

**C21-3900** Safety Burner

## Accessory

**C21-3950** Foot Switch




C21-2700 & C21-3000



C21-3200



C21-3300

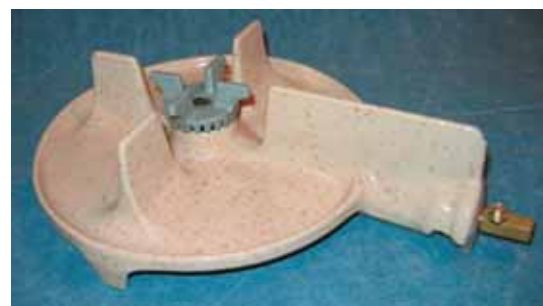

C21-3400



C21-3500



C21-3610



C21-3900